


Features & Benefits

- Anterior-to-posterior bow fits the Pediatric Femur, providing a stable, anatomically appropriate construct for mid-shaft fractures
- All plates will accommodate a non-locking or fixed angle locking screw in every hole, allowing for compression and providing a stable, locking construct
- Broad selection of sizes accommodates a wider range of patients
- Provides a superior alternative to Spica Casting
- Allows for percutaneous insertion
- Provides increased stability over Casting or Flexible Nails
- Facilitates quicker healing and return to function
- Simple easy to use instrumentation

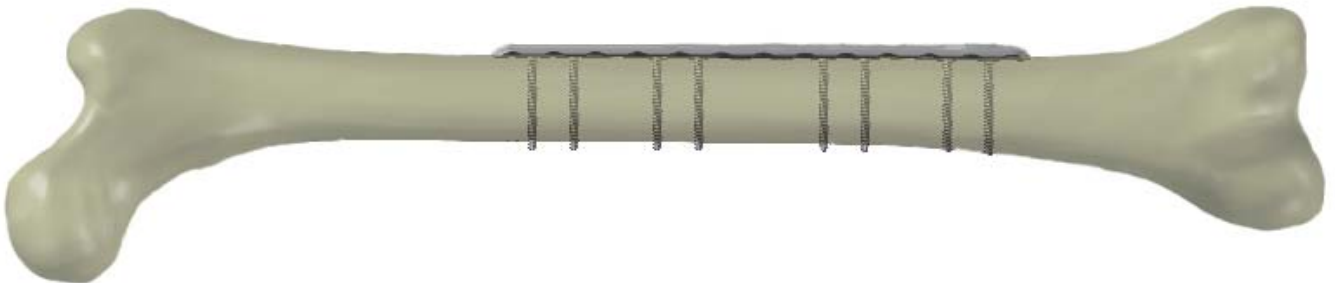

3.5 Plates

- Sizes available: 10, 12, 14, 16 & 18 Holes
- Plates accept 3.5 mm locking cortical screws and 4.0 mm non-locking cancellous screws in full and partial threads

4.5 Plates

- Sizes available: 10, 12, 14, 16 & 18 Holes
- Plates accept 4.5 mm locking cortical screws




Features & Benefits

- Metaphyseal Contour and anterior-to-posterior lateral bow fits the Pediatric Femur and eliminates the need for intraoperative contouring
- Metaphyseal flare and fixed angle screws allow for stable fixation in tough distal femur fractures without endangering the growth plate
- All plates will accommodate a non-locking or fixed angle locking screw in every hole, allowing for compression and providing a stable, locking construct
- Broad selection of sizes accommodates a wider range of patients
- Provides a superior alternative to Spica Casting
- Allows for percutaneous insertion
- Provides increased stability over Casting or Flexible Nails
- Facilitates quicker healing and return to function
- Simple easy to use instrumentation

- Plate is contoured distally to fit the metaphysis
- Metaphyseal flare allows for bi-cortical fixation in multiple planes distally

•Screws are fixed parallel to the axis of the growth plate

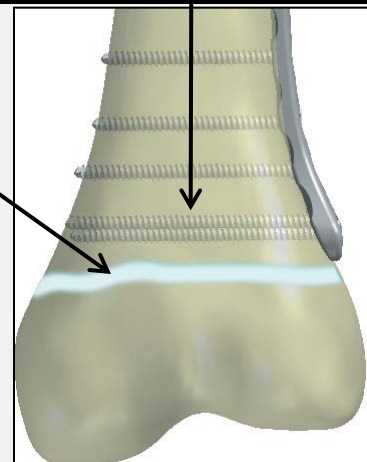

3.5 Plates

- Sizes available: 10, 12, 14, 16 & 18 Holes
- Plates accept 3.5 mm locking cortical screws and 4.0 mm non-locking cancellous screws in full and partial threads

4.5 Plates

- Sizes available: 10, 12, 14, 16, 18 & 20 Holes
- Plates accept 4.5 mm locking cortical screws